

Data sheet

Drying Oven SLW 115 Smart



The photo above is for reference only, may show additional options not included in standard equipment. The real appearance, particularly color and structure of the material may differ from the ones presented in the photo.

Advantages of the SMART controller:

- 4,3", clear, full colour touch screen
- LAN, USB ports for data transfer
- multi-segment time and temperature programs
- visual and sound alarm
- internal memory for programs and data storage
- event registry
- user manual for direct download
- Quick change of program parameters
- Alarm Bar
- operating with gloves on



Smart - preview screen






A. POLOK – KOWALSKA, S. KOWALSKI 44-300 WODZISŁAW ŚLĄSKI, UL. KOKOSZYCKA 172 C

TECHNICAL DATA	
TECHNICAL DATA	
air convection	forced
chamber capacity [I]	112
working capacity [I]	112
controller	microprocessor PID
display	4,3" full colour touch screen
TEMPERATURE	
temperature range [°C]	5°C above ambient temperature +300°C
temperature resolution every [°C]	0,1
temperature fluctuation at 105°C [+/-°C]*	0,2
temperature variation at 105°C [+/-°C]*	2,0
over temperature protection	class 2.0 to DIN 12880 / class 3.1 (option)
CHAMBER	
door type	solid / door with viewing window (option)
interior	
Smart	DIN 1.4301
IG Smart	DIN 1.4301
housing	
Smart	powder coated sheet
IG Smart	stainless steel linen finish

overall dims	[mm] /1/
--------------	----------

Overall alling [mini] /1/	
width A	660
height B	850
depth C	710
internal dims [mm]	
width D	460
height E	540
depth F	450

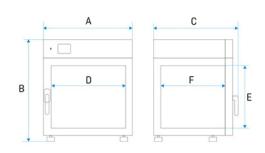

shelves (standard max)	2 7
max shelf workload [kg] /2/	25
- reinforced shelf version (PW) [kg] /3/	50
max unit workload [kg]	60
- reinforced unit version (W) [kg] /4/	120
weight [kg]	65



A. POLOK – KOWALSKA, S. KOWALSKI 44-300 WODZISŁAW ŚLĄSKI, UL. KOKOSZYCKA 172 C

ELECTRICAL PARAMETERS	
voltage**	230V 50-60Hz
nominal power [W]	2500
warranty	24 months
manufacturer	POL-EKO-APARATURA

all the above technical data refer to standard units (without optional accessories)

* - fluctuation measured in centre of the chamber; in space, variation (K) calculated for chamber as:

K = +/- (T average max. - T average min.) / 2

- ** other power supplies on request
- 1 depth doesn't include 50 mm of power cable
- 2 on uniformly loaded surface
- 3 reinforced shelf
- 4 reinforced version

OPTIONS AND ACCESSORIES

	Order number: */A	door with viewing window
INOX	Order number: */P INOX	Stainless steel wire shelf INOX
	Order number: */PP	Perforated shelf
	Order number: */PW	Reinforced shelf
	Order number: KUW GN*/*	Stainless steel cuvettes
Ø	Order number: QLK*	Castors
	Order number: */S or */S INOX	Table with castors
	Order number: */ST or */S INOX	Base on castors
KG	Order number: */W	Reinforced version
	Order number: OWW/OWW LED	Interior lighting



Order number: HEPA

HEPA

Bank: ING BŚ o/Wodzisław Śl. Nr rachunku: 90 1050 1403 1000 0022 0879 5498 Nr rej. W BDO 000005609

HEPA Clean Air Filter

NIP: 647-216-26-87 REGON: 276111984 Sąd Rejonowy w Gliwicach, X Wydział Gospodarczy KRS 0000083696



A. POLOK – KOWALSKA, S. KOWALSKI 44-300 WODZISŁAW ŚLĄSKI, UL. KOKOSZYCKA 172 C

LAB	
DESK	

Order numberLabDesk LabDesk software

>0<

Order number: BRT/*/L or IQ/OQ/PQ Calibration and IQ, OQ, PQ qualification

DIN 3.1

Order number: */3.1 Over temperature protection 3.1 class according to DIN 12880



Order number: BPP 12 Battery backup for display

ALARM 🔲 🗖

Order number: PORT ALARM Dry alarm contact



Order number: USBK USB cable



Order number: LANK LAN cable



Order number: OCZ/20 Non-standard access port 20 mm



Order number: OCZ/30 Additional access port 30 mm



Order number: OCZ/60 Non-standard access port 60 mm



Order number: OCZ/100 Non-standard access port 100 mm



Order number: KD Access control

